

Be Sure of Your Steam Quality

Steam Separator/Filter Kit



Application

The STERIS Steam Separator/Filter can be used with all steam sterilizers and steam generators. The steam separator/filter is installed in the steam supply line for each piece of equipment.

Description

Steam Separator/Filters remove solid and liquid contaminants which may be present in the steam supply line. Maximum steam capacity of 300 lb./hr (can be parallel-piped). The removal efficiency is 1.3 micron in steam service. Utilizes 1" NPT steam line connection.


Benefits

- Helps eliminate condensation in steam lines that may contribute to "wet packs".
- Removes solid and liquid contaminants.
- Helps remove contaminants that contribute to spotting of loads.
- Supplied gauge allows for monitoring of steam pressure.
- Available in "standard" length or "double" length.

Reference List/Order Numbers

Steam Separator/Filter

Standard Length - 168 lb./hr @ 50 PSI	(DIM. A = 13 1/8")	P-760265-859
Double Length - 300 lb./hr @ 50 PSI	(DIM. A = 23 1/4")	P-760265-860
Kit, Element - 10"		P-760265-861
Kit, Element - 20"		P-760265-861

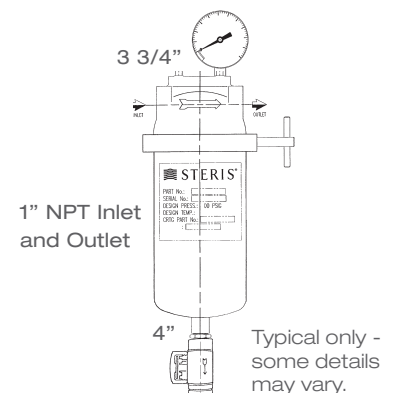

NOTE: Because of STERIS's continuing program of research and development, all specifications and descriptions are subject to change without notice. Obtain certified drawings for design and installation.

This data is intended for the exclusive use of STERIS customers, including architects or designers. Reproduction in whole or in part by others is prohibited.



STERIS Corporation
5960 Heisley Road
Mentor, OH 44060-1834 · USA
440-354-2600 · 800-548-4873
www.steris.com