

Septihol™ Ethyl WFI st Alcohol Solution kills a range of Gram-positive and Gram-negative bacteria, and has fungicidal activity against Trychophyton interdigitale. Its EPA biocidal and fungicidal claims are proven in presence of 5% organic soil on hard, nonporous surfaces. The formula contains 70% (v/v) reagent-grade alcohol and 30% (v/v) United States Pharmacopeia (USP) Water For Injection (WFI), with no stabilizers. Septihol Ethyl WFI st Alcohol Solution products are double bagged and gamma irradiated for easy transfer into aseptic areas.



Septihol™ Ethyl WFI st Alcohol Solution

Ideal for microbial control, residue removal and glove decontamination

Features

EPA biocidal and fungicidal claims on hard, nonporous surfaces in the presence of 5% organic soil

Ready to use without any mixing

Formulation does not contain stabilizers

Double bagged, gamma irradiated by a validated process

0.2 micron filtered, filled in Cleanroom Grade B (Class 100)

Each lot USP sterility tested

Benefits

Facilitates use as part of a qualified microbial control program

Easy to use and saves the time and labor of disinfectant preparation

Eliminates errors in dilution preparation and the need to validate preparation procedures

Evaporates quickly with least amount of residue


Allows for easy entry into controlled environments

Minimal amount of particulates for use in controlled environments

Certificates of Analysis contain irradiation and sterility results

Applications

- Glove sanitization
- Surface disinfection and microbial control of:
 - » Countertops
 - » BSL Hoods
 - » Carts and Other Equipment

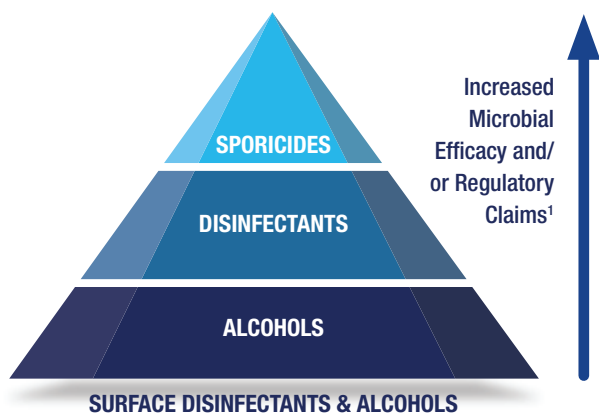

¹ Products that fall into the categories at the bottom of the pyramid are most frequently used and are generally not sporicidal. Progression up the pyramid indicates stronger performance overall and a broader spectrum of claims.

An important aspect to any microbial control program is including rinsing on some frequency to control residues. Septihol Ethyl WFI st Alcohol Solution is ideal for this important application. Its biocidal and fungicidal activity provides an additional level of microbial control during the final step of disinfecting aseptic areas. Formulated without surfactants or stabilizers, the use of Septihol Ethyl WFI st Alcohol Solution as part of a rinsing routine will maintain a clean, residue-free environment. Our Technical Service group can help customize your contamination control program to take advantage of this excellent product.

Septihol Ethyl WFI st Alcohol Solution Ordering Information (United States)

Description	Quantity	Part Number
22 oz (650 mL) trigger sprayer	12 per case	1S1177
32 oz (946 mL) trigger sprayer	12 per case	1S1196
1 gallon (3.78 L) bottle	4 per case	1S1108