



Reliance™ 100 Series

Laboratory Glassware Washers

Reliance 100 Series Washers are high-quality systems specifically designed for research intensive laboratory settings. They offer automated cleaning and drying of typical laboratory glassware, as well as miscellaneous plastic items used in research facilities.

Unmatched Flexibility

The Reliance 100 Series includes three different models to address a diverse range of laboratory research cleaning challenges.

RELIANCE 100

Primarily designed to hold general glassware such as beakers, utensils, test tubes, petri dishes and other wide-mouth items.

RELIANCE 100LS

Specially designed to clean and dry narrow neck glassware such as volumetric and Erlenmeyer flasks and pipettes.

RELIANCE 100XLS

Designed for cleaning and drying of narrow neck glassware in contamination sensitive or critical applications.



Release latches allow upper racks to be easily removed for unloading or exchanging with another rack.

More Flexibility

All three models of the Reliance 100 Series are available in under counter and freestanding styles with an optional 17-inch (43.2-cm) high base stand.



The new color touch screen operator interface allows operators to select and store up to 200 factory-set or user-set programs. The control also provides remote monitoring capabilities and cycle data can be exported via USB/Ethernet ports.



The Reliance 100LS includes a lower spindle rack for efficient cleaning of narrow neck glassware with options for an upper rack.




The side compartment of the Reliance 100XLS houses a 99.97% efficient HEPA filter, automatic liquid detergent and rinse aid dispensers.




The Reliance 100XLS has advanced features for controlling and monitoring contaminants, including water conductivity and HEPA-filtered drying.

Maximum Cleaning Efficacy


The Reliance 100 includes upper and lower standard racks to accommodate a variety of accessory baskets and inserts to ensure efficient cleaning and drying of most typical laboratory glassware.

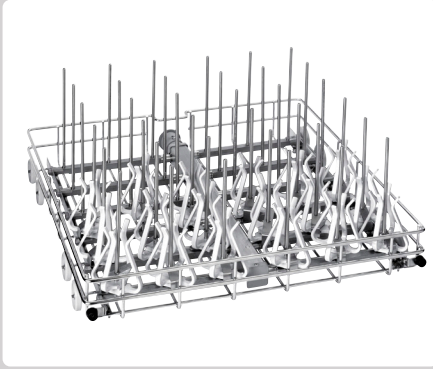

A 36-spindle rack is included in every Reliance 100LS and XLS model. It can also be used in the Reliance 100 version.

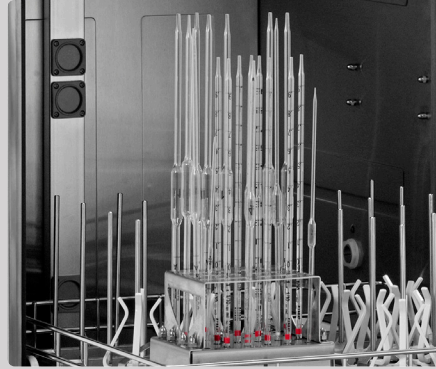

Pipettes may be thoroughly cleaned and dried in Reliance 100LS and XLS Washers.

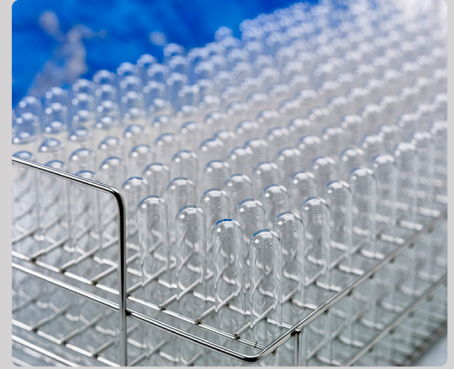

Several inserts are available to ensure complete solution coverage for various types of glassware.



Support

Research solutions from STERIS encompass far more than equipment alone. Every system is backed by service support from an unmatched direct force of technicians worldwide. From specialized chemistries through facility planning, installation and ongoing operational support, STERIS stays with Customers throughout their entire research cycle to maximize uptime and operational life of their equipment.